



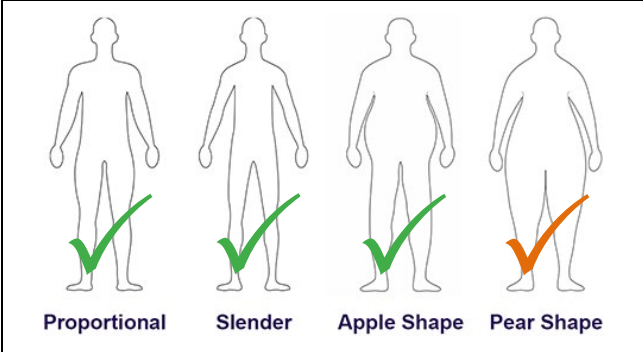


**CAUTION:** This document is intended as a guide only and in no way replaces an individual patient risk assessment. Patient specific conditions, functional independence measure (FIM) and weight distribution may influence the most suitable sling size as well as the most suitable type.

 <p>Proportional    Slender    Apple Shape    Pear Shape</p>	<p>Body shape and size are important considerations when prescribing a sling. This quick reference tool is designed to help determine the compatibility of an individual with this type of sling.</p> <p>✓ = Good general compatibility.          ✓ = May be suitable.          ✗ = Probably not compatible. <b>Use only with extreme caution.</b></p>
---	--

### Key Dimension: Top of Head > Seat

Take the client measurement as shown in the diagram below and reference it against the corresponding measurements in the colour-coded chart to help determine the most suitable sling size.

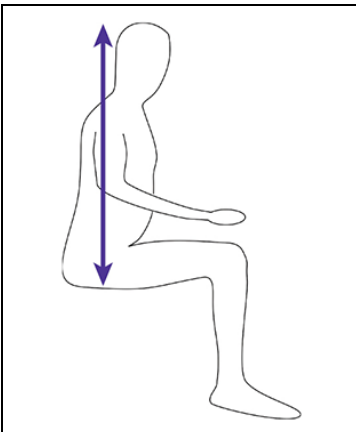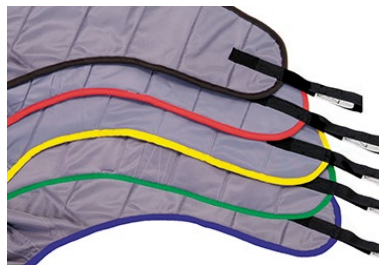	Sling Size	Dimensional Range		<p>Oxford slings have colour-coded binding matching this chart to help quickly identify the size.</p> 
		Min	Max	
	XS	610mm	720mm	
	S	700mm	810mm	
	M	790mm	900mm	
	L	880mm	990mm	
	XL	970mm	1080mm	

For further assistance, core dimensions for each sling size are provided on the next page.

Selecting a sling for an individual involves more than simply measuring their size. Other considerations, as part of a comprehensive risk assessment, should include, but are not limited to:


- Type of transfer (e.g. bed to chair, chair to toilet, etc)
- Level of support required
- Patient conditions
- Physical and cognitive ability of the patient
- Abilities and limitations of the carer
- Compatibility of the sling with other equipment
- Is a commode aperture required?
- Is the sling to be used for bathing?
- Potential obstacles and hazards during transfer



**It is critical to assess all possible risks when planning to lift and transfer a patient from one surface to another. Failure to assess the risk correctly may result in serious injury to the patient and/or caregiver.**


## Core Sling Dimensions: Oxford UltraFine Deluxe Leg

UltraFine Deluxe Leg			
Sling Size	Dimension A	Dimension B	Dimension C
XS	790mm	530mm	860mm
S	880mm	610mm	970mm
M	970mm	650mm	1040mm
L	1070mm	710mm	1140mm
XL	1180mm	800mm	1250mm

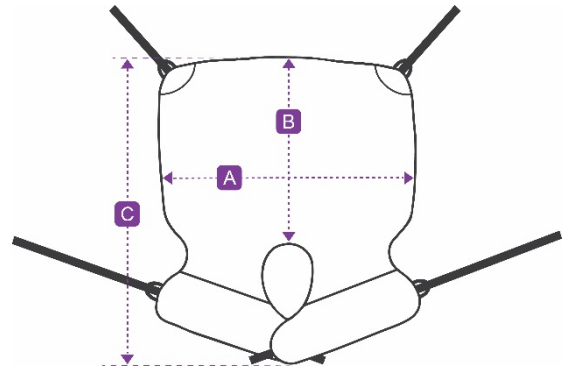

\*All dimensions are subject to manufacturing tolerances.

UltraFine Deluxe Leg with Head Support			
Sling Size	Dimension A	Dimension B	Dimension C
XS	790mm	780mm	1100mm
S	880mm	890mm	1240mm
M	970mm	980mm	1370mm
L	1070mm	1070mm	1470mm
XL	1180mm	1170mm	1620mm

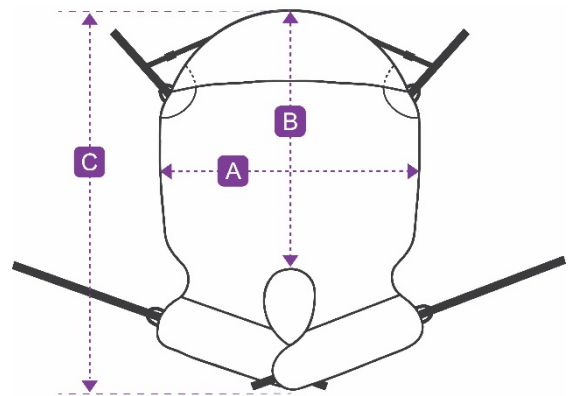

\*All dimensions are subject to manufacturing tolerances.

### Custom Sling Service:

Joerns Healthcare offer a custom sling service should you have any non-standard or patient specific requirements. Please contact Joerns Healthcare for more information and guidance. Extended lead times apply.