



Domus 1

A cost effective solution to reducing the risk of pressure ulcers

APEX
Care for a Healthy Life

Domus 1 is ideal for prevention and treatment of patients at low risk of pressure ulcers. This lightweight and compact alternating pressure pad is composed of 130 6.3 cm high interlocking bubble cells, and can be used as an alternative to heavier products such as gel overlay mattresses.

The Domus 1 is an overlay system designed to be used in conjunction with a full thickness foam mattress, it is simply installed on top of the mattress to increase the levels of support for the patient. It works by the individual bubbles of air inflating and deflating to help reduce the amount of pressure being exerted on the patients skin at any one time.

The system is ideal for short term use and easy to install, the Domus 1 comes with an easy to use manual pressure adjustable pump and securing extension flaps, fitted to the head and foot end of the bubble pad to prevent the pad from sliding. The Domus 1 is designed to be a great first line defence against pressure ulcers in both the acute and community sectors.




Simple intuitive pump design

Mattress Specification	
Mattress Dimensions	: 196x90x6.4 cm
Mattress Weight	: 2.3 KG
Securing flaps at head and foot end	
6.3 cm Bubble Pad Overlay	

Pump Specification	
Maximum Patient Weight	: 100 KG
Simple intuitive design	
Manual weight adjustment dial	
Pump Size	: 23.7x11.2x9.5 cm
Pump Weight	: 1.4 KG
Bubble Cycle Time	: 10 minutes
Voltage	: AC230V, 50/60Hz, 120V, 60Hz
ABS Cover	: ABS Flame Retardant Pump Cover
Hanging Brackets	
Extra Long Power Cable	

Pressure Ulcer Classification Guide

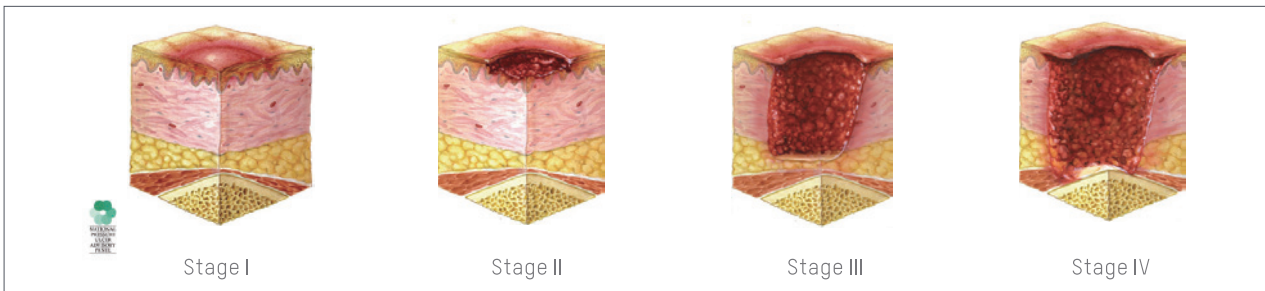

HEADQUARTERS

Apex Medical Limited

Unit 21, Optima Park,
Thames Road, Crayford,
Kent, DA1 4QX
T: 01322 520 560
F: 01322 554 182
enquiry@apexmedicalcorp.co.uk



www.apexmedicalcorp.co.uk