



Propad Revolve

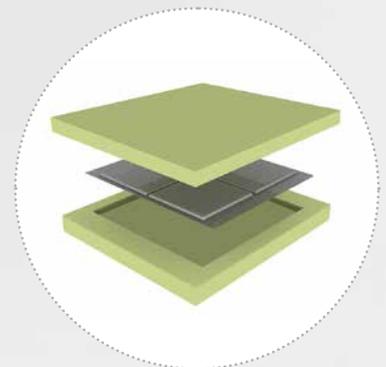

Flat, high risk (V) and very high risk (Si)
cushion that can be placed in any direction

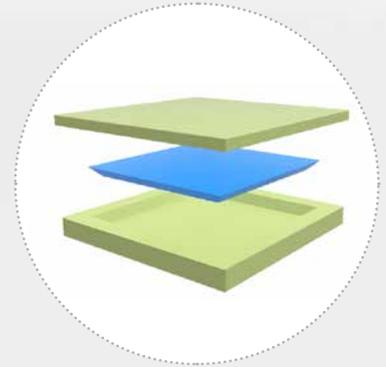

Introducing the **Propad Revolve** cushions

- ▶ Flat cushion
- ▶ Up to 30 stone (190kg)¹
- ▶ Suitable for high to very-high risk
- ▶ Visco and Silicone options
- ▶ Sewn or welded wipedown covers
- ▶ Omnidirectional



High Resilient moulded Foam

- ▶ The **Propad Revolve** cushions are manufactured from the finest quality materials. The same high resilient foam is used within the clinically respected **Matrx Flo-tech** range.




Omnidirectional

- ▶ Versatile and efficient with a design that allows for omnidirectional use from the top, bottom, front or back. Offering increased flexibility and saving time on seating preparation.




Cover options

- ▶ Sewn: blue and black
- ▶ Welded: blue

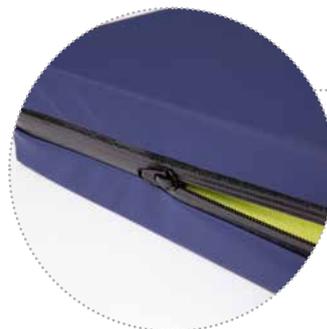

Intelligent fabric

- ▶ A sewn cover comes as standard which features a concealed zip and a full length flap. Completely welded seams available as an option for maximum ingress protection.

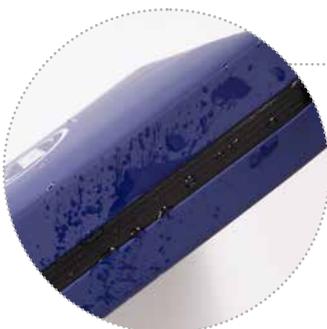


Aquaguard zip

- ▶ An innovative new waterproof repellent zip is included with the welded seamed cover option. Unlike standard zipped covers, the superior Aquaguard withstands an extensive shower test and does not absorb moisture.

¹. Referring to the caveat on page 5 referencing the size of cushion.

Invacare **Propad Revolve Si** (very high risk)

4 EXCLUSIVE FEATURES

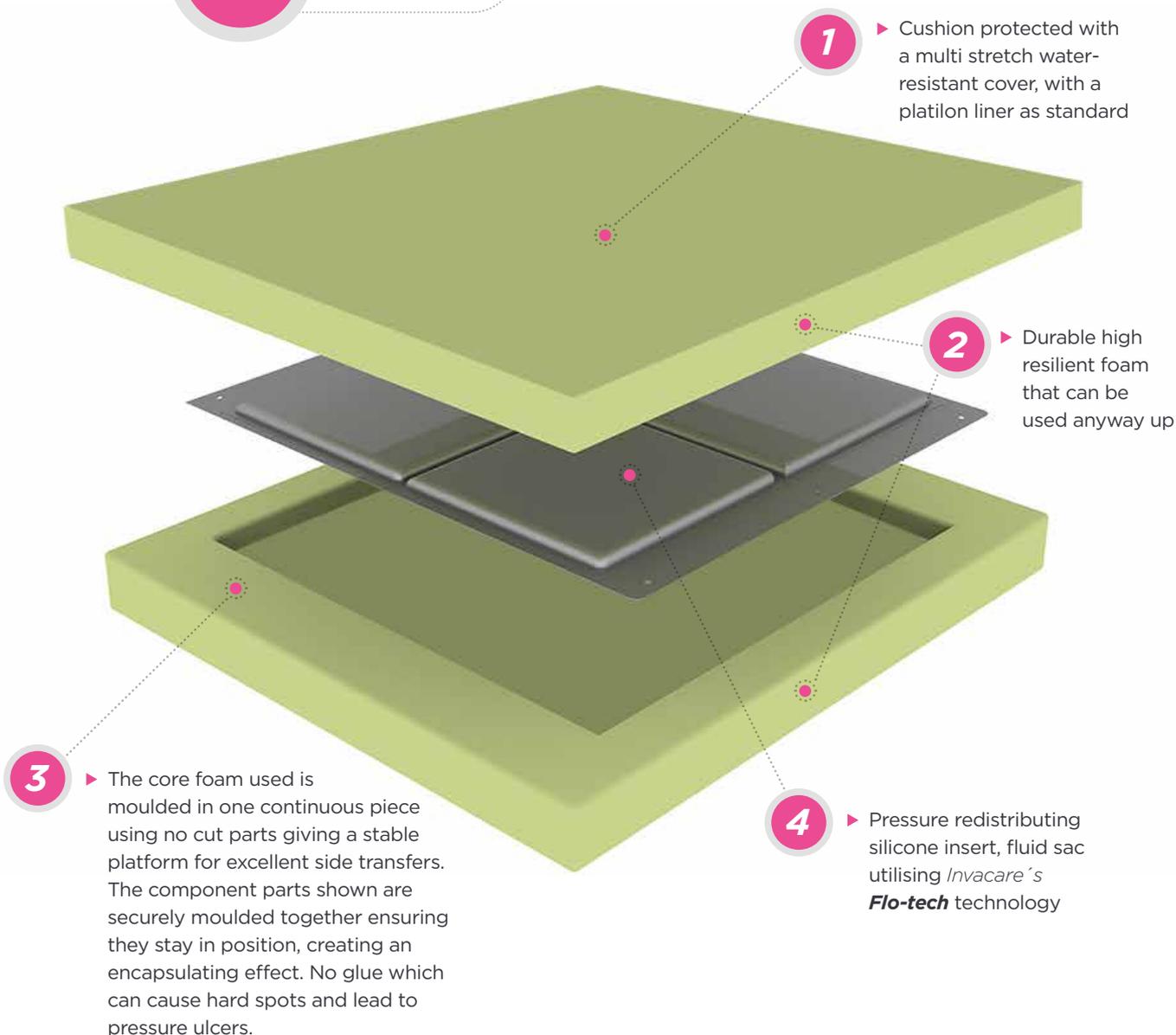

Additional features

▶ Fluid and foam advantages

The embedded, non-temperature sensitive silicone (Si) fluid sac envelops protruding bony prominences for superior skin protection and shear management.

▶ High resilient foam

Top quality moulded HR-foam provides high-level immersion for excellent pressure redistribution.

Invacare **Propad Revolve V** (high risk)

5

**EXCLUSIVE
FEATURES**

1

▶ Lightweight comfortable foam that gently moulds to patient shape, providing excellent comfort and pressure redistributing properties

2

▶ Visco Elastic insert that provides excellent immersion

3

▶ Durable high resilient foam that can be used anyway up

4

▶ Cushion protected with a multi stretch water-resistant cover, with an optional platilon liner

5

▶ One-piece moulded design is excellent for side transfers

Additional features

▶ **Visco-elastic inner layer**

The **Propad Revolve V** combines an outer-layer of durable moulded high resilient foam with a viscoelastic foam inner layer in a lightweight package.

▶ **Design features**

The angled base design ensures seamless transition while the visco layer provides high level immersion and pressure redistribution and stabilization.

Technical data ▼

For more information about this product, including the product's user manual, please visit www.invacare.co.uk



Height ▼



Total product weight ▼



Max. user weight ▼

	Height ▼	Total product weight ▼	Max. user weight ▼
Propad Revolve Si	7,5 cm	2.75 kg	36cm to 41cm - 100kg 43cm to 51cm - 125kg up to 159kg 56cm to 61cm - 190kg
Propad Revolve V	7,5 cm	2 kg	36cm to 41cm - 100kg 43cm to 51cm - 125kg to 159kg 56cm to 61cm - 190kg



Max. cleaning temperature ▼



Materials used ▼



Warranty ▼

	Max. cleaning temperature ▼	Materials used ▼	Warranty ▼
Propad Revolve Si	80°C	Foam, Fluid	3 years
Propad Revolve V	80°C	Foam, Visco	3 years

1. Maximum user weight limits are set provided that the cushion selected is of appropriate width and depth for the intended user/s.

Cushion Sizes ▼

		Width						
		36	38	41	43	46	48	51
Depth	cm							
	inch	14	15	16	17	18	19	20
	36	■	■	■				
	38	■	■	■				
	41	■	■	■		■		
	43				■	■	■	
	46			■		■		■
51							■	